



## High-Power Xenon Light Source Installation and Operation Manual



For Product: HPX-2000-HP-DUV

Document: 000-10000-081-02-201503



**WARNING**

---

Protective Eye Wear Must Be Worn When  
Using This Instrument -  
Intense Ultraviolet Radiation Present  
See Important Safety Notices inside.

---

A HALMA COMPANY

## AMERICAS & WORLD HEADQUARTERS

---

**Phone:** +1 727-733-2447

**Fax:** +1 727-733-3962

**Sales:** info@oceanoptics.com

**Orders:** orders@oceanoptics.com

**Support:** techsupport@oceanoptics.com

### **Ocean Optics, Inc.**

830 Douglas Ave.  
Dunedin, FL 34698  
USA

### **Manufacturing & Logistics**

4301 Metric Dr.  
Winter Park, FL 32792  
USA

## EUROPE, MIDDLE EAST & AFRICA

---

**Phone:** +31 26-319-0500

**Fax:** +31 26-319-0505

**Email:** info@oceanoptics.eu

Germany : +49 711-341696-0

UK : +44 1865-811118

France : +33 442-386-588

### **Sales & Support**

Geograaf 24  
6921 EW Duiven  
The Netherlands

### **Manufacturing & Logistics**

Maybachstrasse 11  
73760 Ostfildern  
Germany

## ASIA

---

**Phone:** +86 21-6295-6600

**Fax:** +86 21-6295-6708

**Email:** asiasales@oceanoptics.com

Japan & Korea: +82 10-8514-3797

### **Ocean Optics Asia**

666 Gubei Road  
Kirin Tower Suite 601B  
Changning District  
Shanghai  
PRC, 200336

---

[www.oceanoptics.com](http://www.oceanoptics.com)

---

### **Copyright © 2012 Ocean Optics, Inc.**

All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, by any means, electronic, mechanical, photocopying, recording, or otherwise, without written permission from Ocean Optics, Inc.

### **Trademarks**

All products and services herein are the trademarks, service marks, registered trademarks or registered service marks of their respective owners

### **Limit of Liability**

Every effort has been made to make this manual as complete and as accurate as possible, but no warranty or fitness is implied. The information provided is on an "as is" basis. Ocean Optics, Inc. shall have neither liability nor responsibility to any person or entity with respect to any loss or damages arising from the information contained in this manual.

---

# Important Safety Notices

1. Do not remove or modify any installed safety device on this equipment. Doing so will void your warranty and create an unsafe operating environment.
2. Dangerous voltages are present in this device. There are NO user serviceable parts inside.
3. Only allow qualified personnel to service this unit.
4. Do not use the unit if it is damaged in any way. Contact your dealer for repair or replacement information.
5. Always screw in the fiber optic cables before starting the instrument.



**WARNING**

Protective eyewear **must** be worn when using this equipment. Intense ultraviolet radiation is present.  
**Never look directly into the light beam**, as this can cause eye damage.

---

## Warranty

Ocean Optics warrants to the original user of this instrument that it shall be free of any defects resulting from faulty manufacture of this instrument for a period of 12 months from the original date of shipment. There are no warranties for the Xenon Bulb (REPAIR-HPX-1) and the Xenon-Bulb-Module (HPX-2000-HP-DUV-BM).

**This instrument should not be used for any Clinical or Diagnostic purposes.** Data generated in these areas is not warranted in any way by Ocean Optics Germany GmbH. Any defects covered by this Warranty shall be corrected either by repair or by replacement, as determined by Ocean Optics Germany GmbH.

There are no warranties that extend beyond the description herein.

This Warranty is in lieu of, and excludes, any and all other warranties or representations expressed, implied, or statutory, including merchantability and fitness, as well as any and all other obligations or liabilities of Ocean Optics including, but not limited to, special or consequential damages. No person, firm, or corporation is authorized to assume for Ocean Optics. Any additional obligation or liability not expressed provided for herein except in writing duly executed by an officer of Ocean Optics:

Ocean Optics Germany GmbH  
Maybachstraße 11  
D-73760 Ostfildern

Tel.: +49 (0)711 34 16 96-0 • Fax.: +49 (0)711 34 16 96-85  
e-mail: info@oceanoptics.de



---

# Table of Contents

About This Manual .....	iii
Document Purpose and Intended Audience.....	iii
Document Summary.....	iii
Product-Related Documentation .....	iii
Upgrades .....	iv
<b>Chapter 1: Setup.....</b>	<b>1</b>
Overview.....	1
Unpacking the HPX-2000-HP-DUV .....	2
Contents .....	2
Connecting the Fiber Optic Cable.....	2
<b>Chapter 2: HPX-2000-HP-DUV Specifications .....</b>	<b>3</b>
Operating Environment .....	3
HPX-2000-HP-DUV Components.....	3
Front Panel .....	4
Rear Panel.....	5
Specifications .....	6
Parts List.....	6
Pinout Information .....	7
Pinout Diagram.....	7
<b>Chapter 3: Operating Instructions.....</b>	<b>9</b>
Operating the Xenon Lamp.....	9
Starting Main Power .....	9
Starting the Lamp .....	9
Turning the Lamp Off.....	10
Warming Up the Lamp.....	10
Operating the TTL Shutter .....	10
Automatic Operation.....	10

Manual Operation .....	11
Operating the Filter Slit .....	11
<b>Chapter 4: Troubleshooting .....</b>	<b>13</b>
<b>Appendix A: Bulb Replacement .....</b>	<b>15</b>
Replacing the Bulb .....	15
<b>Index .....</b>	<b>17</b>

---

# About This Manual

## Document Purpose and Intended Audience

This document provides you with an installation section to get your system up and running.

## What's New in this Document

This version of the *HPX-2000-HP-DUV High-Power Xenon Light Source Installation and Operation Manual* updates the part numbers.

## Document Summary

Chapter	Description
Chapter 1: <a href="#">Setup</a>	Contains a list of package contents and unpacking instructions. Also contains instructions for connecting the fiber optic cable.
Chapter 2: <a href="#">HPX-2000-HP-DUV Specifications</a>	Contains operating environment specifications, as well as other physical details of the product, a parts list, and pinout information for the 15 DB-15 connector.
Chapter 3: <a href="#">Operating Instructions</a>	Provides instructions for operating the Xenon lamp and the TTL shutter.
Chapter 4: <a href="#">Troubleshooting</a>	Contains a table of troubleshooting information.
Appendix A: <a href="#">Bulb Replacement</a>	Provides instructions for changing the bulb.

## Product-Related Documentation

You can access documentation for Ocean Optics products by visiting our website at <http://www.oceanoptics.com>. Select *Technical Operating Instructions*, then choose the appropriate document from the available drop-down lists.

# Upgrades

Occasionally, you may find that you need Ocean Optics to make a change or an upgrade to your system. To facilitate these changes, you must first contact Customer Support and obtain a Return Merchandise Authorization (RMA) number. Please contact Ocean Optics for specific instructions when returning a product.



---

# Chapter 1

# Setup

## Overview

The following sections provide instructions on unpacking and setting up your HPX-2000-HP-DUV High-Power Light Source.

Before using the HPX-2000-HP-DUV for the first time check for transport damage. Be sure to adhere to all warnings on the unit and in this manual.



**HPX-2000-HP-DUV High-Power Light Source**

# Unpacking the HPX-2000-HP-DUV

## ► Procedure

1. Unpack your lamp assembly carefully. Although the lamp is rigidly mounted, dropping this instrument can cause permanent damage.
2. Inspect the outside of the instrument and make sure that there is no damage. Do not use the instrument if damage is present. Contact your dealer for repair or replacement information, if necessary.
3. Use this instrument in a clean laboratory environment (see [Operating Environment](#)).

## Contents


Your HPX-2000-HP-DUV package should contain the following:

- ❑ HPX-2000-HP-DUV High-Power Light Source
- ❑ Power cord
- ❑ UV safety goggles
- ❑ One IC-DB15-2 interface cable for shutter operation
- ❑ This manual

## Connecting the Fiber Optic Cable

### ► Procedure

To connect the fiber optic cable to the HPX-2000-HP-DUV,

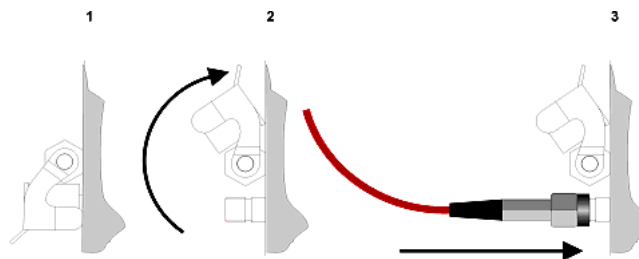

1. Locate the cap on the front of the HPX-2000-HP-DUV.
2. Lift the cap on the front of the HPX-2000-HP-DUV to expose the SMA connector.
3. Connect the fiber optic cable to the SMA connector.

---

## Chapter 2

# HPX-2000-HP-DUV Specifications

This section provides information on the operating environment, physical controls, and dimensions of the HPX-2000-HP-DUV. It also provides pinout information for the connector.

## Operating Environment

The following table provides information on optimizing the operating environment of your HPX-2000-HP-DUV.

Operating Environment	The HPX-2000-HP-DUV Unit . . .
Moisture	Is designed for operation in dry rooms only.
Ventilation	Should be situated so that its location or position does not interfere with proper ventilation.
Heat	Should be situated away from any device that emits excessive heat.
Object and Liquid Entry	Should be positioned so that objects do not fall on top of the unit. Additionally, ensure that no liquids are spilled into the enclosure through openings.

## HPX-2000-HP-DUV Components


The following sections describe the components located on the front and rear of the HPX-2000-HP-DUV unit.

# Front Panel



Component	Description
Power LED	Displays state of main power to lamp.
Xenon ON/OFF	Press and hold the ON/OFF switch for more than 1 second to turn the lamp on or off.
Status LED	Status of the Xenon lamp: <ul style="list-style-type: none"> <li>• Lit – Xenon lamp is on</li> <li>• Blinking – Lamp is in setup mode or cool down mode. The Lamp cannot be started when the Status LED is blinking.</li> <li>• Off – Xenon lamp is off</li> </ul>
Shutter Switch	OPEN - Shutter open CLOSE - Shutter closed TTL - operation by external TTL- signal (HIGH=OPEN, LOW=CLOSE)
Mechanical Protection – SMA Connector	Covered to protect users from unintentionally looking directly at the beam of light. Used for SMA connections only. Connect the fiber cable to the HPX-2000-HP-DUV BEFORE turning the lamp on to avoid unnecessary exposure to UV radiation. Avoid direct contact with the fiber end.  Always wear proper eye protection when using the HPX-2000-HP-DUV lamp.
Four Recessed Screws	Unscrew to remove the unit's housing. See <a href="#">Bulb Replacement</a> .

# Rear Panel

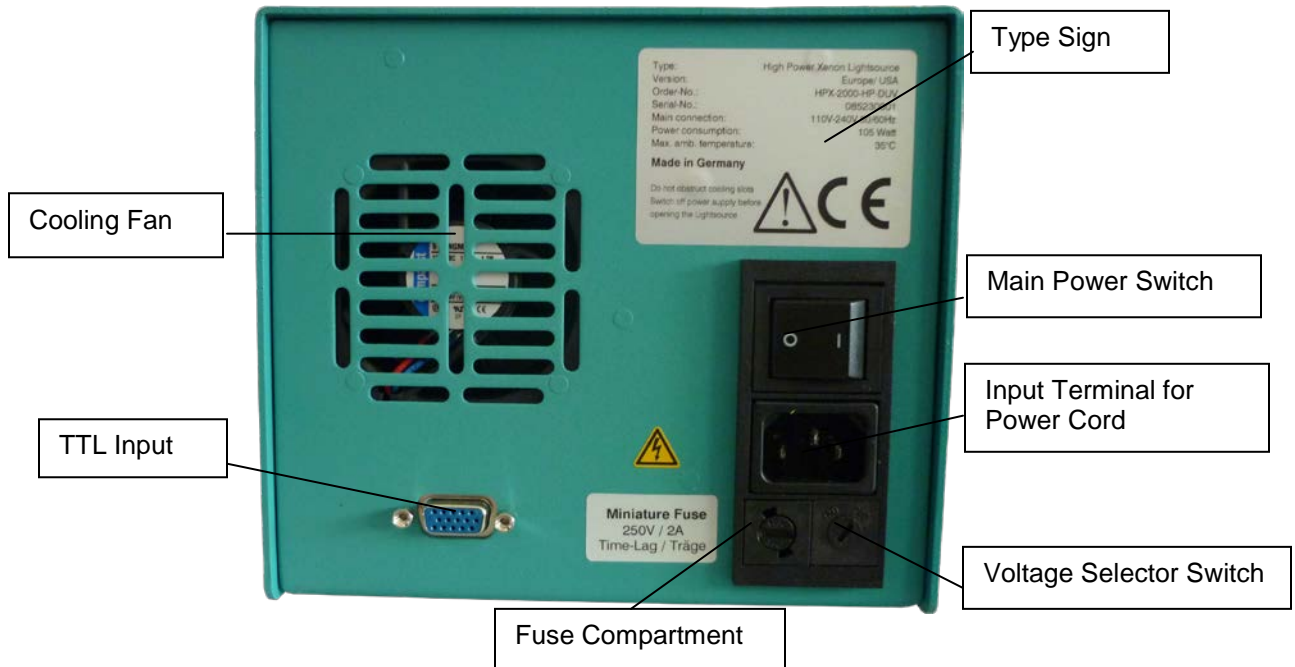

Component	Description
Voltage Selector Switch	Select regional voltage setting for lamp (AC 230-240V, 50/60Hz or AC 110-115V, 50/60Hz).
Main Power Switch	Turn on to supply power to the HPX-2000-HP-DUV. The Power LED lights when this switch is in the ON position.
Fuse Compartment	Contains the fuse to protect the unit against overload: Fuse Type: Miniature fuse 5 x 20 mm, 2.5 Amp slow blow
Input Terminal for Power Cord	Plug power cord into this terminal. <b>Note:</b> Only connect the power cable to the lamp when the Main Power Switch is in the OFF position.
TTL Input	SUB-D 15-pin connector for automatic shutter control.
Type Sign	Label providing lamp information: <ul style="list-style-type: none"> <li>- Type High Power Xenon Lightsource</li> <li>- Version Europe / USA</li> <li>- Order-No. HPX-2000-HP-DUV</li> <li>- Serial-No. Xxxxxx</li> <li>- Mains connection 115 / 230V 47-63Hz</li> <li>- Power consumption 115 Watt</li> <li>- Max. Ambient Temperature 35°C</li> </ul>
Cooling Fan	Cools the interior of the HPX-2000-HP-DUV. Do not obstruct.

# Specifications

Specifications	Criteria
Spectral Range	185 – 2000 nm
Light Output Stability	~ 1%
Warm-up time	25 - 40 minutes
Performance guaranteed temperature	5°C – 35°C
Humidity	5 - 95% without condensation at 40°C
Power Consumption	105 W
Current	5.9 A/DC maximum
Optical Power in 600 $\mu$ Fiber	> 150 mW
Optical Fiber Connection	SMA (optional FC/PC; SC...)
Lifetime Average / Guaranteed	2000h / 1000h
Shutter-Input	TTL-Input, up to maximum 5 Hz
PIN position at SUB-D-15 pin	Shutter PIN 13: TTL / PIN 10: Ground
Input Line	115 / 230VAC
Markings, Directives	CE; VDI/VDE 0160; EN 61010
Lamp Input Current	5.4 A
Weight	5kg
Dimensions	150 mm x 170 mm x 260 mm

## Parts List

Description	Part Number
High Power CW Xenon Light Source with Shutter	HPX-2000-HP-DUV
Replacement Lightsource Module complete with bulb housing optics and shutter, for plug and play replacement on customer side	HPX-2000-HP-DUV-BM
Trade in of used HPX-2000-HP-DUV-BM for Xenon bulb exchange at the supplier side	REPAIR-HPX-1


Modification of device specifications and designs to improve performance may occur without notice.

# Pinout Information

The following table contains pinout information for the HPX-2000-HP-DUV Light Source:

Pin	Description
1	na
2	na
3	na
4	na
5	na
6	na
7	na
8	na
9	na
10	Ground
11	na
12	na
13	TTL Signal – Shutter control
14	na
15	na
na = not applicable	

## Pinout Diagram

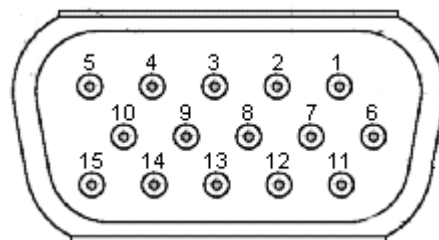





# Operating Instructions

## Operating the Xenon Lamp

The following sections provide instructions on operating the Xenon lamp in the HPX-2000-HP-DUV Light Source. The HPX-2000-HP-DUV unit must be in a horizontal position for it to work.

### Starting Main Power

► **Procedure**

To apply main power to the unit,

1. Ensure that the proper voltage is set on the Voltage Selector Switch on the rear panel of the unit (see [Rear Panel](#)).
2. Ensure that the Main Power Switch is OFF. Then, connect the power cord to the Terminal Input on the rear panel of the unit.
3. Turn the Main Power Switch ON. The Power LED lights, indicating that the unit is receiving power. The Xenon bulb remains off after power-on.

### Starting the Lamp

Hold the Xenon ON/OFF switch down for 1 full second to illuminate the HPX-2000-HP-DUV lamp.

The lamp cannot be started when the Status LED is blinking. The Status LED blinks for about 5 seconds after the main power is switched on (main power switch is located on rear of unit).



**Protective eyewear must be worn when using this equipment - Intense ultraviolet radiation present. Never look directly into the light beam, as this can cause eye damage.**

---

## Turning the Lamp Off


Turn the Xenon lamp off by holding the Xenon ON/OFF switch down for one full second. The Status LED blinks for about 30 seconds while the lamp cools down.

## Warming Up the Lamp

The HPX-2000-HP-DUV requires 25-30 minutes of operation to reach a state of thermal equilibrium. During this warm-up period, the intensity of the output power can vary substantially.

If applications require extreme intensity stability, the lamp should be warmed up for an additional 30-45 minutes. Once warmed up for this amount of time, the lamp will reach specified drift values.

The following graph shows a typical behavior of the intensity of a HPX-2000-HP-DUV during warm-up.

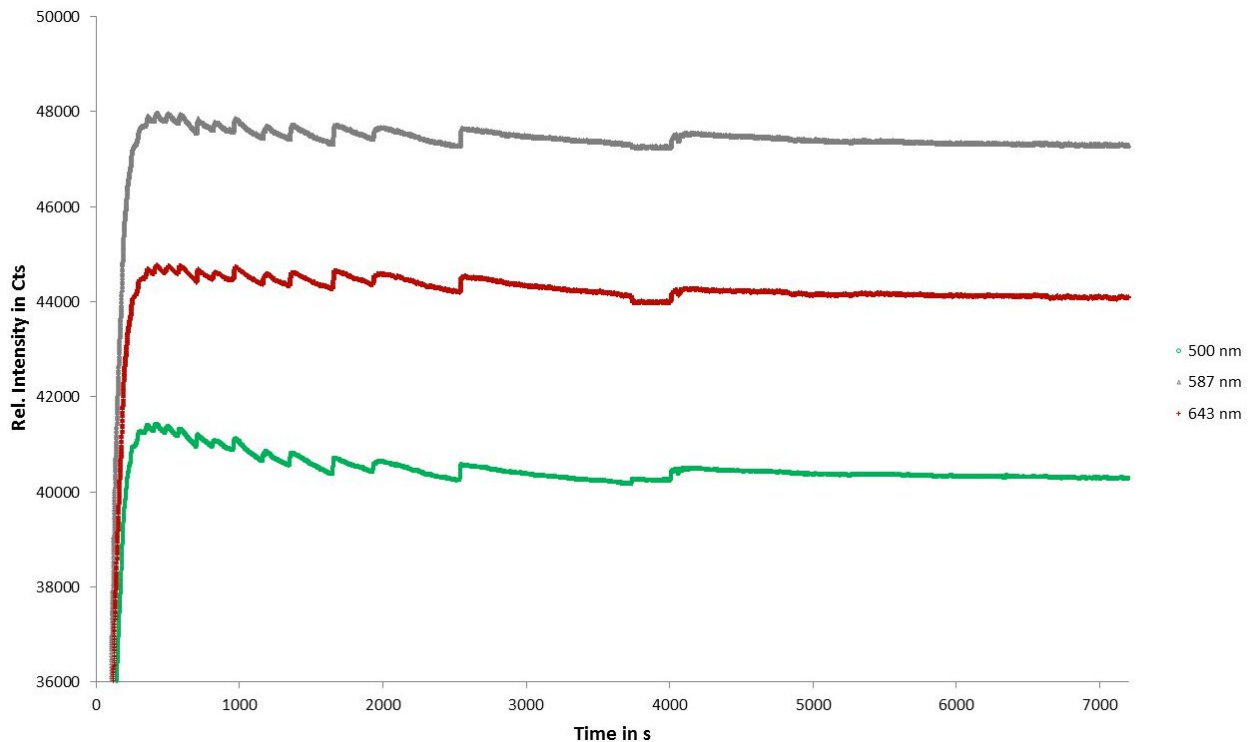

## Operating the TTL Shutter

### Automatic Operation

For automatic operation, plug the SUB-D 15-pin TTL connector into the appropriate socket on your spectrometer.

## Manual Operation

Manually set the operating mode of the HPX-2000-HP-DUV with the Shutter Switch as follows:

- OPEN - Shutter open
- CLOSE - Shutter closed
- TTL - Controlled over external TTL signal (HIGH=OPEN, LOW=CLOSE)

## Operating the Filter Slit

### ► *Procedure*

1. Rotate the light beam protection cap to open the filter slit (see [Front Panel](#)).
2. Insert your filter with a maximum size of 1” round or square into the filter slit.
3. Rotate the light beam protection cap to close the filter slit.



# Troubleshooting

If the power supply or lamp does not seem to functioning properly, check the following:

Issue	Probable Cause	Resolution
Power switches on, but no LEDs light.	Line power not present	<ol style="list-style-type: none"><li>1. Check line voltage</li><li>2. Check if the lamp is plugged in properly.</li></ol>
	Fuse defective	Check fuse
Xenon lamp does not start	Lifetime of the lamp is exhausted	Replace the lamp module
	Lamp's internal connection plug is not closed properly	Open unit (see <a href="#">Bulb Replacement</a> ) and close connector plug.



---

## Appendix A

# Bulb Replacement

To order replacement bulbs for the HPX-2000-HP-DUV, order item number HPX-2000-HP-DUV-BM.

## Replacing the Bulb

---

### WARNING

**During operation, the lamp reaches a temperature of approximately 250°C. Before changing the lamp's bulb, please wait at least 20 minutes for cooling before touching the lamp module.**

---

#### ► *Procedure*

Follow the steps below to change the bulb in the HPX-2000-HP-DUV:

1. Turn the Main Power Switch OFF.

---

### Caution

**Connection/disconnection of the power cord from the Xenon lamp MUST only be done when the Main Power Switch is turned OFF. Because the output of the power supply is not galvanically separated from the line voltage, the connecting lead for the Xenon lamp can carry a voltage >42 V when the Main Power Switch is turned ON.**

---

2. Disconnect the power cord from the unit.
3. Remove the 3 stub screws (Figure 1, #2) using the 1.5 mm Allen wrench included with the replacement bulb.
4. Remove the cap (Figure 2, #3) and the beam cover (Figure 2, #4) from the HPX-2000-HP-DUV.

**A: Bulb Replacement**

5. Remove the four recessed screws (Figure 2, #1), both front and back on both sides of the HPX-2000-HP-DUV.
6. Remove the cover housing of the HPX-2000-HP-DUV.



2


**Figure 1**



4

**Figure 2**

7. Loosen the screws (Figure 3, #5) by hand, disconnect the plug (Figure 3, #6), and remove the old bulb from the HPX-2000-HP-DUV.
8. Place the new HPX-2000-HP-DUV-BM bulb inside the housing, connect the plug, and replace the screws.
9. Close the housing and screw the 4 screws (Figure 2, #1) into the sides of the housing.
10. Replace the cap (Figure 2, # 3) and the beam-cover (Figure 2, #4).

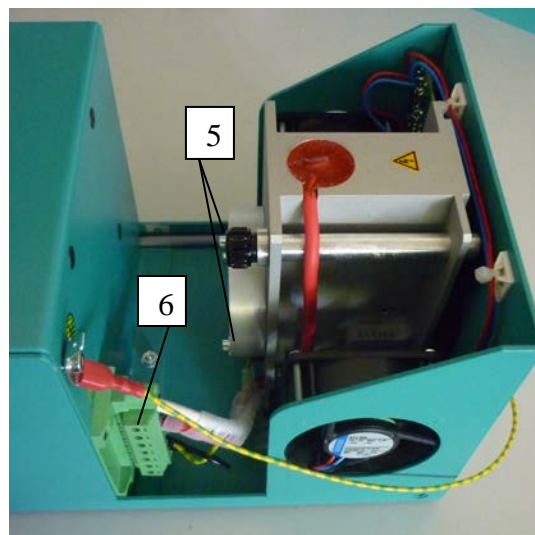

**Figure 3**



---

# Index

## B

bulb replacement, 15

## C

cable connection

    fiber optic, 2

components, 3

    front panel, 4

    rear panel, 5

connecting the fiber optic cable, 2

## D

document

    audience, iii

    purpose, iii

    summary, iii

## F

fiber optic cable

    connecting, 2

filter slit, 11

## O

operating environment, 3

## P

package contents, 2

parts list, 6

pinouts, 7

product-related documentation, iii

## S

setup, 1

specifications, 3, 6

## T

troubleshooting table, 13

TTL shutter, 10

    automatic operation, 10

    manual operation, 11

## U

unpacking procedure, 2

upgrades, iv

## W

warranty, A

## X

Xenon Lamp, 9

    starting, 9

    turning off, 10

    warming up, 10

